

MA.X2 MOISTURE ANALYZERS



release date 03-02-2017




MA.X2 moisture analyzers are hi-tech measuring instruments intended for quick determination of relative moisture content, dry mass content and other parameters in samples of different substances.

Highlights of MA X2 series moisture analyzer:


- ż Modern interface
- ż **DRYING FORECAST – Prognosis of the Drying Process Result**
- ż Touchscreen
- ż Programmable display
- ż Easy configuration of information arrangement
- ż Wide range of databases – work optimization
- ż Databases – security and archiving of your work results
- ż USB interface allowing export of results and databases
- ż E2R software allowing on-line control over your work

Interface


Hot-Keys

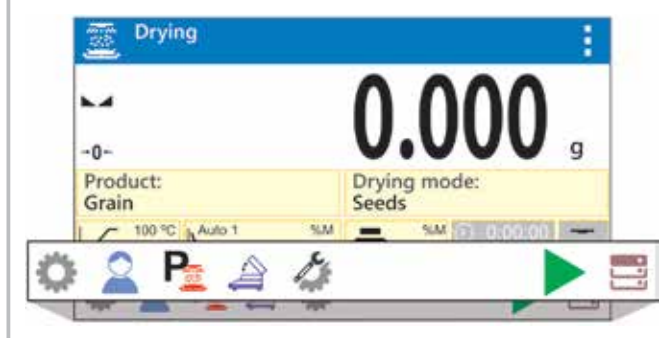

Labels and information fields

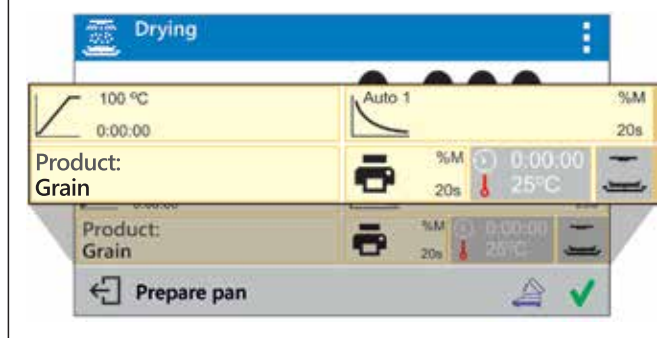

Functions keys

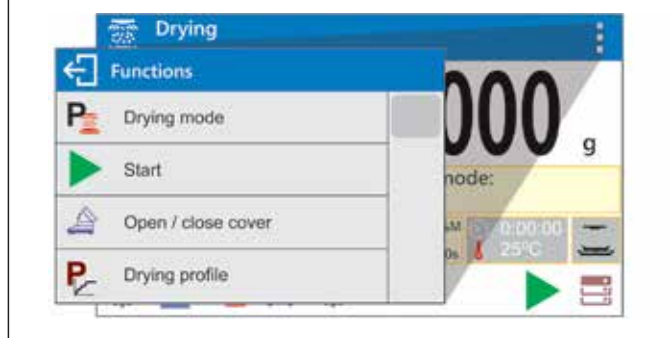

Databases

Drying process can be carried out for any temperature and set for a specific product using following databases:

PRODUCTS	5000	USERS	100
DRYING PROGRAMS	200	DRYING PROCESS REPORTS	5000
CUSTOMERS	100	PACKAGINGS	100

GLP/GMP Printouts

An example of a printout

Information on drying parameters and moisture content result can be printed immediately using following elements: header – GLP – footer. Each of those elements can be customized in order to create a printout consistent with the requirements (GLP/GMP).

Drying process report	
HEADER	
Start date	2016.07.02
End date	14:24:02
Operator	John Smith
Product	Grain
Drying mode	Seeds
GLP	
Drying parameters	
Drying profile	Standard
110 °C	
Finish mode	Automatic3
1mg/60s	
Indication	%M
Interval	30 s
Start mass	2,595 g
FOOTER	
0:00:30	? 0,694%M
0:01:00	? 0,694%M
0:01:30	? 0,694%M
0:07:00	? 0,694%M
0:07:23	? 0,694%M
Status	Completed
End date	2016.08.03
End time	0:07:23
Drying time	14:24:02
End mass	2,315 g
Indication	10,790 %M

DRYING FORECAST - Prognosis of the Drying Process Result

MA X2 series moisture analyzers are now equipped with DRYING FORECAST function that enables to shorten the drying process. If you do not care about the highest accuracy and increased measurement error (ranging from 5% to 20% of the end value) you can shorten the drying process up to 6 times. It is an option dedicated for users who require fast and multiple estimation of materials moisture content.

DRYING FORECAST method is a prognosis of the end result, carried out before the drying process is completed.

Based on characteristics of current drying curve, created online, moisture analyzer estimates the end result of the drying process.

It is an approximate result characterized with small error of +/- 10% of the end result for most of the products and +/- 20% for products that require drying process to be longer (such as plastics, gels and materials containing more than 30% of water).


The calculation algorithm adapts to the drying material features and estimates first end result. The value of the prognosis is more accurate when the process takes more time. You can specify the time and accuracy of the process.

Enabling DRYING FORECAST Function in MA.X2


By default, DRYING FORECAST function is disabled.

In order to enable it:


1. Select 'Settings' menu located in the left bottom corner of the screen

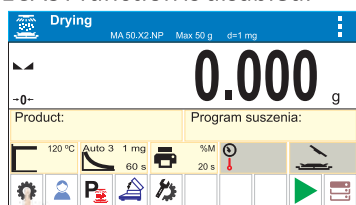

2. Select 'Working modes' submenu

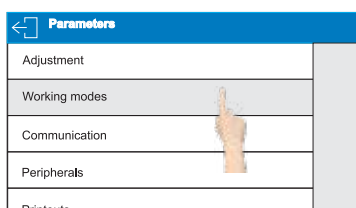


3. Select 'Drying' mode

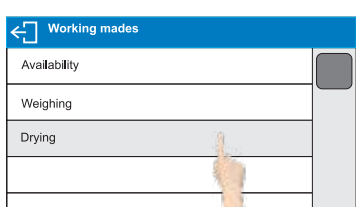


4. Use the scroll bar on the right to reach end of the menu


5. Enable DRYING FORECAST function by selecting 'Prognosis' parameter

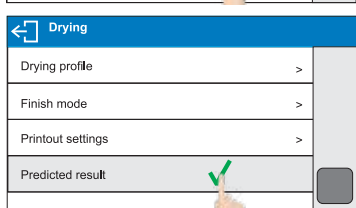

Enabling Prognosis Function

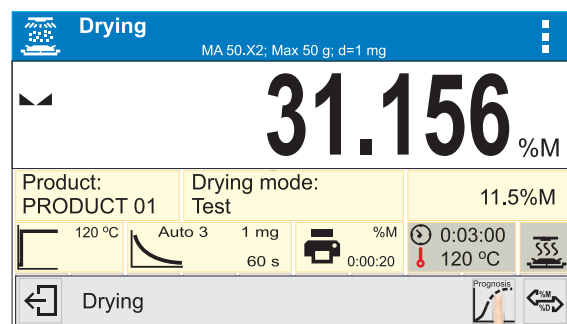

Prognosis of the Drying Process Result

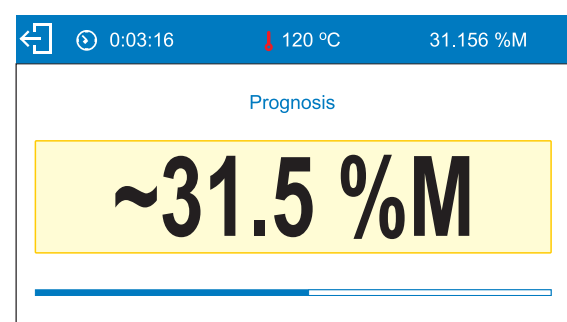

Drying Process Report Containing Prognosis

Drying process report		
Start date	2016.11.14	
Start time	14:41:55	
User	John	
Product	PRODUCT 01	
Drying mode	Test	
Drying process parameters		
Drying profile	Standard	
120°C		
Finish mode	Auto3	
1mg/60s		
Indication	%M	
Interval	20s	
Start mass	0.796 g	

0:00:00	0.000	%M
0:00:20	1.508	%M
0:00:40	7.412	%M
0:01:00	15.578	%M
0:01:20	21.859	%M
0:01:40	26.256	%M
0:02:00	29.020	%M
0:02:20	30.528	%M
0:02:40	31.030	%M
0:03:00	31.156	%M
0:03:16	31.156	%M
Status	Aborted	
Drying time	0:03:16	
End mass	0.548 g	
Indication	31.156 %M	
Predicted result	31.5%M	

Signature


.....

Technical specification:

	MA 50/1.X2	MA 50.X2	MA 110.X2	MA 210.X2
Max capacity	50 g	50 g	110 g	210 g
Readability [d]	0,1 mg	1 mg	1 mg	1 mg
Tare range	-50 g	-50 g	-110 g	-210 g
Max sample weight	50 g	50 g	110 g	210 g
Moisture readout accuracy	0,0001 %	0,001 %	0,001 %	0,001 %
Moisture content repeatability	0,05% (sample weight of 2 g), 0,01% (sample weight of 10 g)			
Max sample height	h= 20 mm			
Weighing pan size	Ø 90 mm, h= 8 mm			
Drying temeprature range	max. 160° C			
Heating module *	IR emitter			
Drying mode	4 drying modes (standard, quick, step, mild)			
Auto switch-off options	4 options (time-defined, automatic, manual, user-defined)			
Additional functions	Control of sample weight prior drying process			
Power supply	230 V			
Display	LCD 5" capacitive touchscreen			
Interface	1×RS 232, 1×USB-A, 1×USB-B, 1×Wireless Module			
Proximity sensors	2			
Net weight/gross weight	4,9 / 6,4 kg			
Packaging size	470×380×336 mm			

* Heating element options: WH - halogen (max= 250° C), NS - metal heater (max=160° C)

Dimensions:

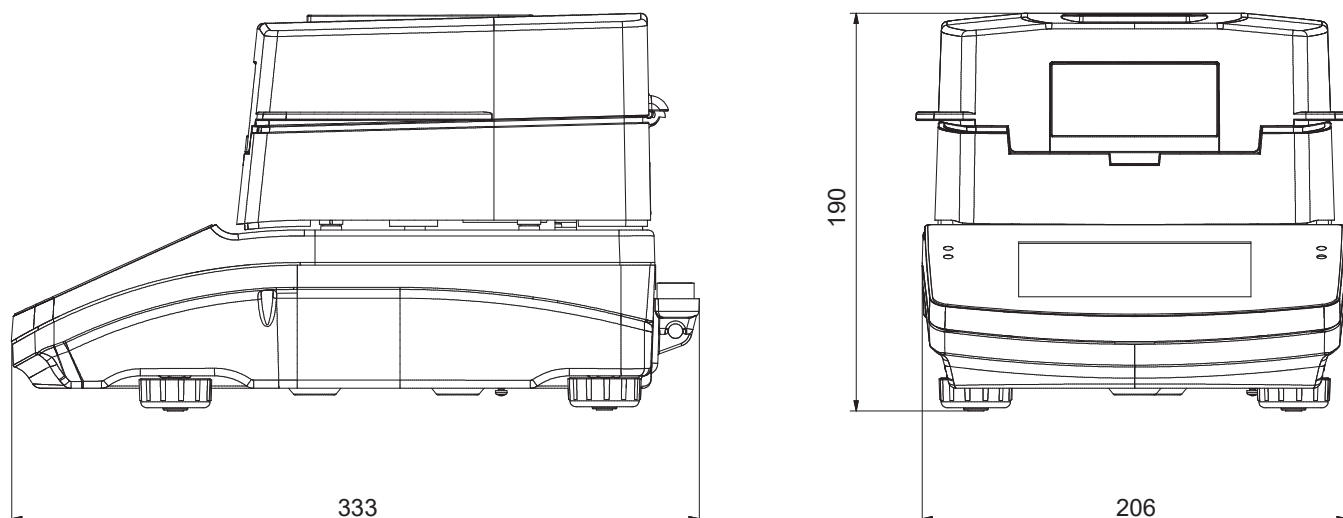

Additional equipment:

Anti-Vibration Table (coated or stainless steel)	GT105k-12/Z Control Thermometer
EPSON Printer	Adjustment Weight
Disposable Weighing Pans	PW-WIN, RAD-KEY PC Software
PC Keyboard	Calbe RS 232 (balance - computer) "P0108"
Water Vapour Permeability Determination Set	Cable RS 232 (balance - Epson, Citizen printer) "P0151"